


MILTEX HI-LIGHT

DOVGAN ENDODONTIC CONDENSERS

NiTi Condensers

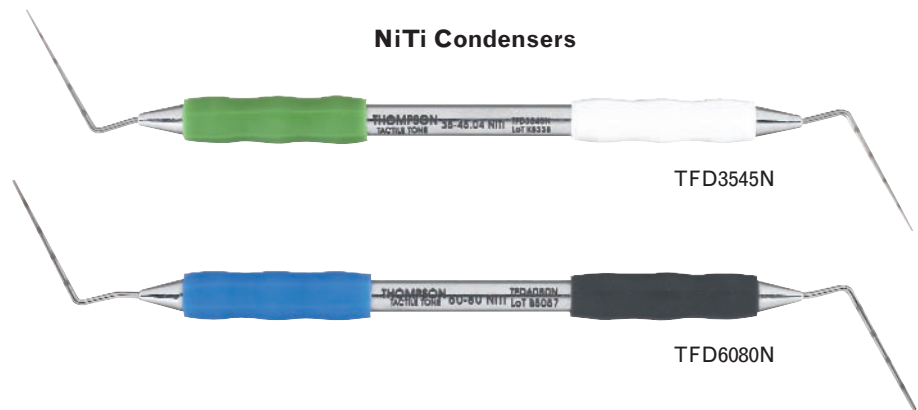

TFD3545N

TFD6080N

Large Condenser



TFD100140X

REF	Description
TFD3545N	DOVGAN Plugger .04 ISO 35-45 Tactile Tone double end NiTi 23 mm (green/white)
TFD3545NXL	DOVGAN Plugger .04 ISO 35-45 Tactile Tone double end NiTi 27 mm (green/white)
TFD6080N	DOVGAN Plugger .04 ISO 60-80 Tactile Tone double end NiTi 23 mm (blue/black)
TFD6080NXL	DOVGAN Plugger .04 ISO 60-80 Tactile Tone double end NiTi 27 mm (blue/black)
TFD100140X	DOVGAN Plugger .02 ISO 100-140 Tactile Tone double end GTX 23 mm (yellow/orange)

The Dovan Endodontic Condensers are specifically designed for vertical condensation. They feature a double-ended design which saves time during procedures and ISO sizes for easier identification.

Sizes 35 – 80 are available in a 23mm version and an 27mm extended length version, and are crafted from NiTi and machined to a .04 taper. This ensures proper compacting of gutta percha around curved canals.

Sizes 100 – 140 are stainless steel with a GTX™ non-stick finish applied to the instrument. The GTX coating helps reduce the sticking and pull back effect gutta percha can have on condensers in the coronal region of the tooth.

All instruments available for individual sale. For additional information, an in-office demonstration, or to place an order, please contact your local Miltex, Inc. authorized Distributor or Miltex, Inc. Customer Service at 1-866-854-8300, email customerservice@miltex.com, or visit us at www.miltex.com.

Miltex
Redefining Excellence

www.miltex.com
toll free 866 854-8300
phone 717 840-9335