


MILTEX HI-LIGHT

DR. ROANE ENDO CONDENSER KIT

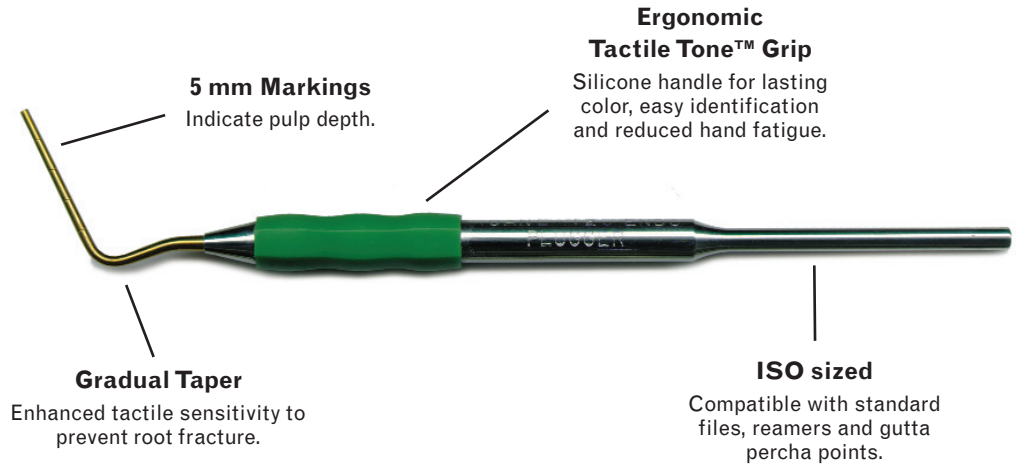


The Dr. Roane Endo Condenser Kit* is specifically designed to improve the handling of gutta purcha and Resilon® during obturation of canals. The smaller sized instruments are made of Nickel Titanium (NiTi) for better compacting of the gutta percha around curved canals. The larger ISO sized condensers have the popular Miltex® GTX™ non-stick coating applied to them. This coating reduces the sticking and pull-back effect gutta percha can have on condensers in the coronal region of the tooth.

The kit includes six condensers from size 45 to 140. All instruments are sterilizable with other stainless steel instruments and are included in a Slimline 7 cassette.

Resilon® is a trademark of Jeneric/Penton, Inc.

*Designed by Dr. James Roane

DR. ROANE ENDO CONDENSER KIT (XFR) CONTAINS ONE OF EACH:

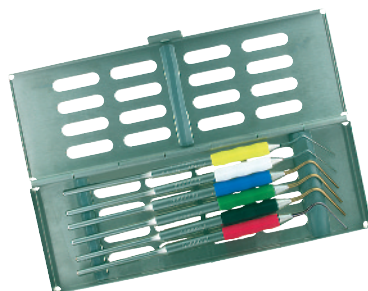

REF	Description
UF45N	45TT SE NiTi Condenser White
UF60N	60TT SE NiTi Condenser Yellow
UF80N	80TT SE NiTi Condenser Red
UF100X	100TT SE GTX Condenser Blue
UF120X	120TT SE GTX Condenser Green
UF140X	140TT SE GTX Condenser Black
3-072007	Slimline 7 Cassette

Miltex
Redefining Excellence

www.miltex.com
toll free 866 854-8300
phone 717 840-9335

For additional information, an in-office demonstration, or to place an order, please contact your local Miltex, Inc. authorized Distributor or Miltex, Inc. Customer Service at 1-866-854-8300, email customerservice@miltex.com, or visit us at www.miltex.com.